

FT Catheter Fixation Tape



Fast and reliable fixation of catheters and tubes

- Designed to secure catheters and tubes to a patient's skin
- Adhesive tape with two resealable wings for anchoring and repositioning of catheters and tubes
- Skin-friendly adhesive provides for reliable and stable fixation for up to 5 days
- Free of Polyvinyl chloride (PVC), Latex, Phthalates/DEHP, Bisphenol and Per- and polyfluoroalkyl substances (PFAS)


CORRECT POSITIONING AND REPOSITIONING

The wings can be opened and closed repeatedly, allowing for repositioning of the catheter or tube.

PATIENT COMFORT

Flexible material conforms to the patient's body, minimising catheter movement and discomfort.

Soft foam reduces pressure points and the risk of pressure injuries.

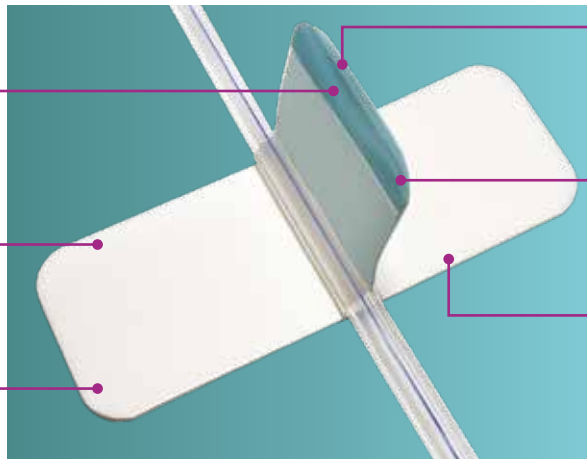

EASY TO USE

Two wings coated with a special adhesive allow for fast securement of catheters and tubes, designed to prevent accidental displacement.

Blue non-stick film with different lengths ensures easy application, even with gloves.

Adhesive PU-foam base coated with a skin-friendly polyacrylate adhesive allows easy and secure application to the patient's skin.

FT Catheter Fixation Tape is suitable for anchoring a variety of devices to the patient's skin, including urinary catheters, nephrostomy tubes, infusion lines, chest drains, central venous lines, wound drains, and peritoneal dialysis tubes. Accommodates catheters and tubes with CH/Fr. 8-36.



Scan to learn more and watch the video



ORDERING INFORMATION

Ref.	Description	Pack QTY.
Q37449	FT Catheter Fixation Tape, 102x38 mm	1/50/600

CONTACT US

UNOQUIP GmbH
Schützengraben 20
8200 Schaffhausen
Switzerland

info@unoquip.com
www.unoquip.com



Scan to enter our website

With reservations for misprints and alterations in product specifications. Printed version not controlled. All trademarks and logos are the property of UNOQUIP GmbH. CE-marked product. © Copyright UNOQUIP GmbH 2025. PM0024 v01MAR2025.